



Wipe Spot Pads


Ultra-absorbent Wipe Spot Pads from Capintec easily absorb any removable contamination for wipe testing. These wipe pads are backed by a card with space for indicating the wipe location, date and time, activity found, and the initials of the user. These soft, highly absorbent cotton pads may be used either wet or dry to discover contamination.

The Wipe Spot Pads are packed 500 to a box. Cases of 5000 are available for facilities that process large numbers of wipes. The adhesive backing allows the wipe to be folded back upon itself and inserted into a test tube for counting.



CAPRAC Printer

An Epson printer is provided as an option for the CAPRAC Wipe Test/Well Counter. This sturdy, high-quality printer provides graphic output with wipe reports, showing the peak energy range of the nuclides on the wipe. The printer utilizes a black/red ribbon to report wipes that are higher than the trigger levels for contamination. Abnormal values are printed in red for easy/fast identification.

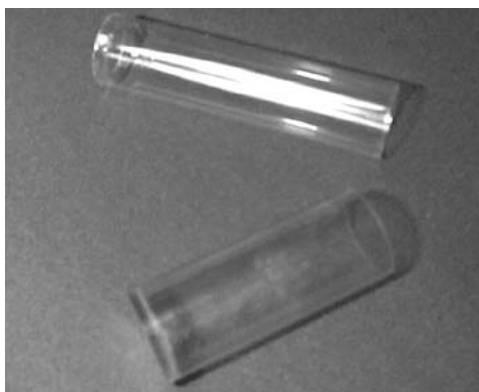

CAPRAC Liner/Test Tube

ITEM #	DESCRIPTION	QTY
5420-0086	Wipes	500
0670-0016	Test Tubes	250
5420-0121	Liners	5
5430-0058	Epson Roll Printer	
5420-0087	Test Tubes	100